

KaVo Services

**DENTAL  
EX  
CE  
LL  
EN  
CE**



The products, equipment and services shown and described in this catalogue are not available in all countries. All information corresponds to the knowledge available at the time of printing.

## Contents

<b>04</b>	<b>Dental Excellence</b>	
	Made in Biberach	04
	Play it safe: use original parts	10
	KaVo Service Partner	14
	KaVo Service Portfolio	16
<b>18</b>	<b>Treatment units</b>	
	KaVo Maintenance	20
	Maintenance kit	22
	Protect+ service package	23
	e-Service	24
<b>28</b>	<b>Treatment unit accessories</b>	
<b>38</b>	<b>Instruments</b>	
	KaVo Original Factory Repairs	40
	KaVo BOX	44
	KaVo BOX PLUS	45
<b>46</b>	<b>Instrument accessories</b>	
<b>50</b>	<b>Material numbers</b>	
	Contact partners	64



# Made in Biberach

At KaVo, dental excellence encompasses not only premium instruments and treatment units, but also excellent service and maintenance protocols. This is the proven level of all-round quality that pays off - every single day.



KaVo has over one hundred years of experience to call upon. It was founded in 1909 in Berlin-Steglitz by Alois Kaltenbach, and since then, all of our products and services have been guided by the motto “Quality and Precision”. We have always set new standards with our innovations, particularly in our core business area of instruments and treatment units. More than 3,000 employees around the world ensure that we provide high quality goods and services, and strive for customer satisfaction at all times. In addition to quality, tradition and sustainability are important attributes at KaVo, we continue to rely on “Made in Germany” quality. Everything comes directly from our main factory in Biberach – from development and assembly of all treatment units to the production and assembly of all instruments, and even our services.







#### Development – made in Biberach

To ensure that new KaVo products meet the strictest requirements, our innovations are always created through teamwork: Specialists from the fields of design, construction, development, marketing and sales collaborate throughout the entire product development process. Initial sketches are converted into 3D data, initial polystyrene models are turned into prototypes, and prototypes eventually into end products.



#### Spare parts production – made in Biberach

Original spare parts are indispensable to preserve the value of the product as well as for safety, and to ensure that your KaVo product remains a genuine KaVo product. Everything – from the manufacturing processes to the production plants – is governed by a single objective: Maximum quality.

#### Instrument repair – made in Biberach

We manufacture and produce everything on-site in Biberach to ensure the safety of the treatment team and the patients, as well as to guarantee the best possible treatment results. KaVo instruments are valuable and their function and quality are intended to last for many years of use. Rather than discarding defective products, it is almost always more lucrative to repair them using original spare parts for KaVo Instruments.

#### Technical support (hotline) and service training – made in Biberach

KaVo-certified specialised dealers always provide high-quality service to our customers. The technical support team is also available to troubleshoot and offer rapid solutions. KaVo Campus Service Training ensures that specialised dealer technicians receive appropriate training. This training not only ensures high KaVo service quality, but also forms an integral component of the required KaVo certification.

#### Customer service – made in Biberach

A high level of emphasis is placed on daily contact with the certified specialised dealer in this regard. Direct and fast communication is the key to ensure that we always provide the customer with what they need. Every day, KaVo Customer Service handles over 500 customer enquiries and orders in German, English, Spanish, French, Italian, Russian, Dutch and Turkish, and ensures that every enquiry is processed within 24 hours.



*Quintessential quality workmanship –  
KaVo instrument sleeves are milled from  
a single piece of raw material.*



# Stay safe – with KaVo original parts

What you should pay attention to when equipment requires servicing.



Medical devices are subject to special laws and standards when compared to many other products that we use in everyday life. These regulations may vary in different countries, but they have one thing in common: They were all created to protect the patient. Product operators or distributors who fail to comply with these regulations – particularly if the patient is endangered or injured as a result – face severe and far-reaching consequences.

More than **50 %**

of instrument failures occur due to inadequate or improper cleaning and care, or due to application errors.



## CONSEQUENCES OF FAILING TO COMPLY WITH THE MANUFACTURER'S SPECIFICATIONS:

- If the shaft is worn, the drill bit holding force might be too low and the drill bit might fall out during treatment
- The shaft could spin in the collet chuck and destroy it
- The part could fall out
- Excessive loads may have been placed on the ball bearings, gears and collet chuck e.g., by tools that are too long
- A grooved/notched shaft may become jammed in the instrument clamping area



**KaVo recommends using only KaVo original spare parts for servicing. Why?**  
KaVo products are not only valuable but supremely robust. They are designed for the user to ensure the best possible treatment quality for many years. A great deal of expertise and care goes into development and production by our KaVo experts. Materials and processes are not only subject to our own demands for the highest quality, but also to strict standards and certifications. The use of KaVo original parts therefore guarantees the same functional quality as a new product, as these correspond exactly to our own strict quality specifications. This ensures that your KaVo product remains a KaVo product.

**KaVo original spare parts are therefore a fundamental requirement in order to preserve the following important benefits. They contribute to...**

- **Durability and value preservation**  
When you use KaVo original spare parts, you will ensure that you are preserving your initial investment, for many years to come. Repairs are done easily and efficiently. This is how we make


an important contribution to the reduction of disposable products, as well as to sustainable and eco-friendly management of increasingly scarce and expensive raw material resources.

- **Functional safety**  
The repaired function must always comply with all KaVo specifications and must not change.
- **Legal safety**  
Medical devices are subject to strict standards and laws. The use of non-original parts contravenes KaVo specifications and can lead to the loss of the standard warranty, or even cancellation of the operating licence. Anyone who sells non-original spare parts, or who uses products repaired with non-original spare parts, or releases such products onto the market should ensure that there is proof of compliance with all laws, standards and certifications. If a patient or a member of the practice team sustains an injury, this will have potentially serious consequences for the product operator and distributor.

By using KaVo original spare parts, you therefore always ensure that you are on the safe side.

# Value preservation of KaVo instruments in the dental practice

Tips for extending the service life of rotary instruments

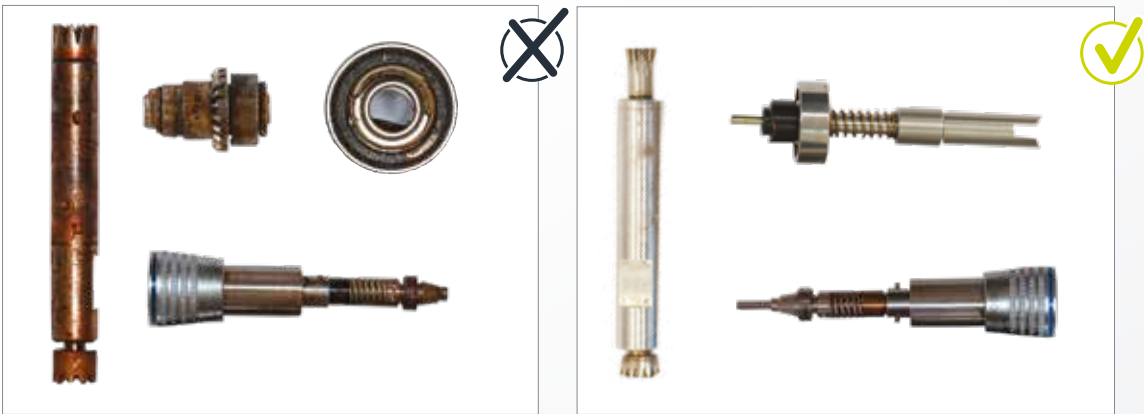

Instruments that have not been correctly cared for   Instruments that have been correctly cared for

## Value preservation of the clamping system

### 1. Rotating instruments/drills

Always ensure that the specifications of the drills that you use have been approved by KaVo. If this is overlooked, the transmission instrument may be damaged. Damage may include defects in the clamping system or the drill shaft.

### 2. Cutter and grinder dimensions

Manufacturer specifications for length, diameter, shaft shape and maximum revolutions per minute are included in the instructions for use for the transmission instrument.

### 3. Removing crowns and bridges

Rotary instrument manufacturers offer crown cutters adapted to different crown materials. For example, they recommend using a special

crown cutter made of hard metal for metal or soft ceramic alloy crowns. In contrast, a diamond crown cutter must be used for zirconium oxide crowns.

Please adhere to the speeds recommended by the drill manufacturer. These speeds are often lower for straight or cross-toothed crown cutters than for diamond crown cutters.

### 4. Relief of the collet chuck

Remove the drills when you store the instruments, since the collet chuck should be relieved during storage to increase its service life. There is also a risk of injury or infection, among others, if instruments are stored with drills.

## Value preservation of ball bearings

### Maintain ball bearings, e.g. with KaVo spray or QUATTROcare:

Inadequate care may result in severe wear and damage of the ball bearings. It is advisable to replace defective ball bearings early to prevent more serious consequential damage.

## Common application errors

### 1. Processing

Never immerse the instrument in a disinfectant/ultrasonic bath because:

- The ball bearings may be destroyed
- Other technical instrument defects may occur

### 2. Instrument care

The service life of the handpieces and contra-angle handpieces or turbines depends greatly on how well the clamping system, ball bearings and the gears are cared for and maintained.

As a rule, automated care using the QUATTRO-care PLUS is preferable to manual care. Please note that KaVo can only guarantee that your KaVo products will work flawlessly if you used KaVo-approved care products. We recommend caring for the instruments after each cleaning and disinfection. We would also advise you to treat the instruments and instrument heads separately.

### 3. Drying and storage

Be sure to dry your instrument thoroughly, both inside and out, to prevent it from being damaged by liquid residues. Remove any liquid residue with compressed air, otherwise corrosion may occur.

## Prevent damage from accidental drops and collisions

As a general rule, take precautions to ensure that your instruments cannot be knocked to the floor easily or dropped, to prevent serious damage. However, if an instrument accidentally falls to floor, you can perform a visual inspection for possible changes or damage to ensure maximum safety.



You can find detailed tips on how to care for your instruments in the "KaVo Care Guide", available at:



[https://p.widencdn.net/bco1zn/BR\\_Pflegesysteme-und-mittel\\_20191017\\_01\\_En](https://p.widencdn.net/bco1zn/BR_Pflegesysteme-und-mittel_20191017_01_En)



# KaVo Service Partner

Proper care, regular maintenance and original spare parts are essential for ensuring that a KaVo product remains a genuine KaVo product. We are here for you – we ensure that we can provide you with excellent service, together with our certified service partners, allowing you to focus all of your attention on your patients during routine work in your dental practice.



1,200

Each year, we train employees from our certified specialised dealers at the KaVo Campus in Biberach. This is how we ensure that your service technician is capable of performing all maintenance and service work efficiently, safely and at the level of quality expected from KaVo. This is to assure you that your investment is optimally protected.

34  
SPECIALISTS

at KaVo support KaVo specialised dental dealers via a technical support hotline, technical documentation and comprehensive training. We take every step to ensure fast and optimal service.

## KaVo Service Partner

### Quality should be close enough to be obvious:

Your KaVo service partner is local to you and is your first port of call for maintaining and servicing your KaVo products.



Genuine parts carry this sign.

KaVo service partners only use KaVo original parts. This is indicated by affixing this symbol.

### Play it safe with KaVo service partners

KaVo service partners receive intensive training from KaVo. They possess these fundamental prerequisites for qualified KaVo service:

- Required training for the technicians
- Required tools and testing equipment
- Exclusive use of KaVo original spare parts

There are so many reasons to perform maintenance and any required repairs of your valuable KaVo products using only KaVo original and spare parts over the long term. It preserves value, as well as user and patient safety. You also save money and protect the environment. On the other hand, the installation of foreign parts may lead to the violation of standards and laws or the loss of warranty, and could decrease the product value, performance and function.



# The KaVo Service Portfolio

KaVo supports its repair partners around the globe with valuable services to ensure that your problems are always solved quickly.



## KaVo Hotline

If you are looking for technical support, this service is available in many countries through the KaVo sales company, or through our certified specialised dealers. For fast and optimal service. The technical hotline at the main plant, located in Biberach-Warthausen, provides additional third-level support at a global level.



## KaVo Campus



Employees and service technicians from our certified specialised dealers are trained at the KaVo Campus by experienced KaVo specialists using the latest training methods.



## Services in the treatment units area



### KaVo Maintenance Kit

For your operational and planning safety.




### KaVo Service Protect+ Package

Extended warranty including other protective extras.



### KaVo e-Service

Installation, data transfer, updates and checks – simple, fast online support for your KaVo equipment.

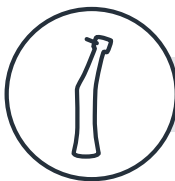

## Service in the instruments area



### KaVo Original Factory Repair\*

You know in advance how much you will save later: Factory repair at attractive and fixed prices. This allows you to know in advance how much you will save later.



### KaVo BOX PLUS\*

Full repair cost protection for 2 years after the repair.

\*This offer is not available in all countries. You can find more information at [www.kavobox.com](http://www.kavobox.com).

# TREATMENT UNITS

THE CENTREPIECE  
OF THE PRACTICE





# KaVo Maintenance: Offering you safety during routine activities at your dental practice

As a medical device manufacturer, we are obliged to stipulate maintenance activities, technical safety inspections and integrated unit tests. This approach serves two purposes: it ensures functional safety and prevents damage, health hazards and legal risks. Regular maintenance also prevents expensive downtimes and helps you to preserve the value of your unit.



KaVo recommends the annual maintenance service for your treatment units to ensure continuous operational readiness, safety and value preservation of the KaVo product. The annual maintenance service also helps to comply with the specifications, standards and laws in your country.

The international standard IEC (DIN EN) 62353 is used to ensure electrical functional safety. It also demands that the functional safety of all technical application and safety-relevant functions are checked and repaired, if required.

Please observe the statutory provisions in your country.

## Dental practice operators are responsible for:

- Conducting the safety check in accordance with IEC 62353
- Complying with hygiene standards

These requirements can be met during the scope of maintenance carried out by KaVo or the certified specialised dealers.

## The following individuals are authorised to carry out repairs and maintenance:

- Technicians from the KaVo branch who have received appropriate product training
- Technicians from KaVo authorised dealers who have been specially trained by KaVo

All certified service technicians are also obliged to only use original KaVo spare parts.



35,000

treatments and

20,000

operating hours will be performed in a treatment unit over an average of

12 years.

During this time, approx.

100,000

litres of water are used.

Source: KaVo empirical values

Maintenance components	Benefits for the dental practice	
Testing of all important functions, e.g. hygiene functions, safety switch, etc.	Functional safety	Safety of the dentists, dental practice team and patients
Cleaning the suction system, Aquamat and selection unit   Replacing wearing parts	Functional safety	Safety of the dentists, dental practice team and patients
Replacing worn parts in the water block   Inspecting for leaks	Functional safety	Safety of the dentists, dental practice team and patients
Inspecting the O-rings on the doctor's and assistant's element, hoses and pressures   Replacing worn parts	Functional safety	Safety of the dentists, dental practice team and patients
Safety inspection in accordance with IEC (DIN EN) 62353	Functional safety	Safety of the dentists, dental practice team and patients

Functional safety  
Safety of the dentists, dental practice team and patients



**Costs time and money:**  
Waiting for service

**Saving you time and money:**  
Planned, periodic maintenance

## The best way to prevent costly downtime: The **KaVo maintenance kit** for your operational safety

The certified technician and the KaVo maintenance kit are an unbeatable team when it comes to ensuring the operational safety of your treatment unit.

The kit supports:

- Safety: Replacement of necessary parts with KaVo original parts
- Efficiency: Systematic workflow of maintenance activities



## Saving you time and money: The extended warranty offered by the **KaVo Service Protect+ package**

The Service Protect+ package offers you the security and support that you need to ensure that you are able to run your dental practice smoothly.

Dental unit and equipment	Safety package and contents
1058/E50 Life with compact or German Technical and Scientific Association for Gas and Water (DVGW) water block	<b>KaVo Service Protect+ 1058/E50 Life</b> <ul style="list-style-type: none"> <li>• Additional two-year warranty <sup>1</sup></li> <li>• 12 L DEKASEPTOL disinfectant and cleaning agent for suction systems <sup>2</sup></li> <li>• 12 L OXYGENAL 6 for dental treatment water germ reduction</li> <li>• 1x footrest protector, to provide simple and aesthetic protection of the upholstery</li> </ul>
E70/E80 Vision with compact or DVGW water block	<b>KaVo Service Protect+ E70/E80 Vision</b> <ul style="list-style-type: none"> <li>• Additional two-year warranty <sup>1</sup></li> <li>• 12 L DEKASEPTOL disinfectant and cleaning agent for suction systems <sup>2</sup></li> <li>• 12 L OXYGENAL 6 for dental treatment water germ reduction <sup>2</sup></li> <li>• 1x footrest protector, to provide simple and aesthetic protection of the upholstery</li> </ul>
1058/E50 Life without DVGW water block	<b>KaVo Service Protect+ comp. 1058/E50Life</b> <ul style="list-style-type: none"> <li>• Additional two-year warranty <sup>1</sup></li> <li>• 6 L DEKASEPTOL disinfectant and cleaning agent for suction systems</li> <li>• 1x complete kit for preventive regular maintenance</li> </ul>

All extended warranties are also available individually, i.e. without the PROTECT+ package. Please see the appendix for more details. The KaVo warranty terms and conditions apply.

<sup>1</sup> After the statutory warranty period of 12 months has expired, free spare parts if required. Consumables and worn parts, as well as labour and travel expenses are not included

<sup>2</sup> Average annual requirement



### Benefits:

- Ensures the legally prescribed levels of hygiene and functional safety
- Reduces downtimes and unplanned costs, which in turn ensures investment and planning security

### Prerequisites

May be purchased up to three months after purchasing a treatment unit The extended warranty General Terms and Conditions can be obtained via your dental dealer. The prescribed annual maintenance is a mandatory condition of purchase.

For more details:

<https://www.kavo.com/accessories-upgrades/protect-service-package-protect>







# Service

Access the best service with just one click.

## KaVo e-Service.

Gain access to four areas with just one click

### 1 e-Installation

Installation, configuration and activation of new software programs

### 2 e-Data

Migration of data to new program versions or a new server

### 3 e-Update

Updating of software programs and installation of extensions with updates

### 4 e-Check

System and error analysis with remote maintenance

## Easy. Consistent. Everywhere. Service when you want it\*

With KaVo e-Service, we offer you simple and fast online support for your KaVo equipment, software programs, applications, databases and server services at attractive and fixed prices.

In many cases, this allows you to spare yourself the effort of scheduling appointments, as well as the technician's travel expenses to and from the site: KaVo IT specialists and service technicians perform the desired services remotely on your dental practice network.

Response times are optimised and downtimes in the practice are minimised because we carry out the KaVo e-Service at a pre-agreed time.

Another benefit is the complete cost transparency and control: Billing is subject to clearly defined and attractive conditions.

The data protection is as secure as the fair conditions: Connection to your systems is established via a specially secured line and therefore offers the highest level of protection. In addition, the connection is only established after you give your express permission.

\*This offer is not available in all countries. Please contact us at [service.multimedia@kavo.com](mailto:service.multimedia@kavo.com)

### Installation of DIAGNOcam Software (single workstation) 1.013.3414

- Checking the system requirements and calculating the duration of the data import from KaVo KiD into the DIAGNOcam software
- Advice on integrating the practice management software, backup and anti-virus protection
- Installation of the DIAGNOcam software on a client and data import from the KaVo KiD software into the DIAGNOcam software.

### Installation of DIAGNOcam on an additional client 1.013.3415

- Checking the client system requirements
- Installation and configuration of the software and connection to the server system

# 46 %

of the service costs can be saved with e-Service, as you are only billed for the effective labour costs. There are also no travel costs, so appointments can be made quickly and downtime is reduced.

### Installation of CONEXIO software (single workstation) 1.013.3410

- Preliminary clarification of the system requirements with the customer or the network administrator
- Discussion of the integration of the dental practice management software, image processing software, backup and anti-virus protection
- Installation of the CONEXIO software on the server, as well as activation and release of the purchased licence
- Configuration of the treatment unit, as well as creation of user accounts and user rights in the CONEXIO software
- Final test of the basic functions and transfer of the pending connection to the practice management and image processing software to the on-site system administrator

### Installation of CONEXIO software, 1x server, 1x client 1.013.3412

- Preliminary clarification of the system requirements with the customers or the network administrator
- Discussion of the integration of the practice management software and image processing software, backup and anti-virus protection
- Installation of the CONEXIO software on the server and on a client, as an example, as well as activation and release of the purchased licences.
- Configuration of the treatment units, creation of user accounts and user rights in the CONEXIO software
- Final test of the basic functions and transfer of the pending connection to the practice management and image processing software to the on-site system administrator

### Installation of CONEXIO 1.013.3413 on an additional client

- Checking the client system requirements as well as installation and configuration of the software and connection to the server system

The demand for  
e-Service performance  
has increased by

# 58 %

in 2020 due to the  
Coronavirus lockdown.

## e-Data

### Moving the database to a new server 1.013.3416

- Analysis of the legacy data to support the server change and calculation of the time needed to transfer the data
- Discussion of the procedure with the customer and/or the network administrator
- Migration of the database to the existing software version and implementation or renewal of licence adjustments, if required (Prerequisite: legacy and new server are fully installed and functionally configured in the same network)

## e-Update

### Update within the main version (server only) 1.013.3418

- Checking the software updates for compatibility with the client systems and, if necessary, informing the customer or network system administrator in preparation for necessary adjustments to devices, licences and operating systems
- Server-side installation of the software updates available free of charge and adaptation based on the previous customer configuration
- Discussion of other update activities to be carried out on-site in the network and instructions on how to do this

## e-Check

### General remote maintenance 1.013.3404

- System or error analysis and technical application support, with billing based on time and effort
- Billing starts once the connection to the connection number provided by the support is established. Billing ends once there is disconnection of the remote maintenance session by the service employee or the customer





# More cost-effective, because it is original: The KaVo Original New Upholstery Set



Even the very best upholstery set will lose its colour intensity or become damaged at some point. If you want to maintain a well-kept and aesthetically pleasing appearance, you can simply renovate your treatment unit instead of replacing it completely. This saves you money and conserves valuable resources. While other replacement upholstery is only refilled with foam and reupholstered, the KaVo original new upholstery set is produced from scratch with the same KaVo original equipment quality, including the upholstery support. You can choose your desired shapes, materials and KaVo colours – and the conditions are so attractive that you won't need to think twice.

*An upgrade to the RELAXline soft cushion is also possible at any time.*



## RELAXline Soft Upholstery: The most beautiful form of softness

Relaxation that can be visualised: Your eyes say: fascinating, elegant colours – graphite, cashmere and agave – combined with a great design. Your senses say: extremely comfortable memo foam, wrapped in extremely soft upholstery material. Your patients will say: I have never been so comfortable, even during longer treatments.



Soft upholstery consisting of memo foam and an imitation leather cover:


- Beautiful, pleasantly soft, leather-like surface
- The upholstery material comes with antimicrobial protection, and is certified in accordance with the United States Environmental Protection Agency (EPA) specifications
- High-quality design seams

- Exclusively handcrafted in a manufactory – Made in Germany
- Available in three colours: Graphite, cashmere and agave, exclusively for the KaVo RELAXline soft upholstery
- The Progress and Comfort backrests for the seat cushion, and double-jointed headrests with mechanical/pneumatic clamping and the motorised headrest






## Upholstery colours

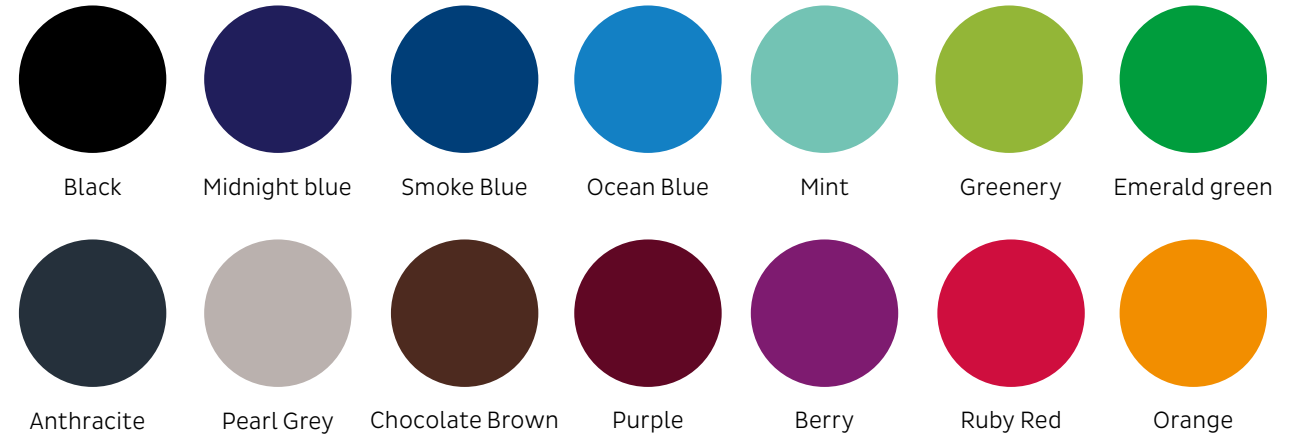

## RELAXline Soft Upholstery

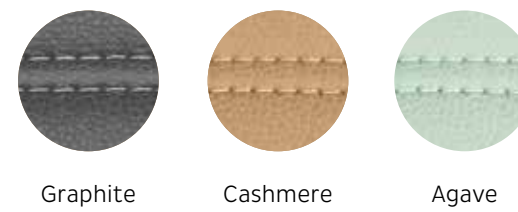

## KaVo Comfort Head Cushion

### Product info

- Comfortable patient positioning, especially during long periods of treatment
- Rotating cushion for maxillary or mandibular treatment
- Easy affixing of the magnetic wedge attachment
- Headrest can be mounted separately: No additional assembly costs
- Heat-sensitive upholstery material for optimal adjustment

### Your benefits

- Quick conversion: with just a flick of the wrist
- Memory foam for comfortable positioning and precise adaptation to any head shape
- Suitable for prophylaxis, regular and long-term treatment



- Cushion easy to turn, depending on whether you are performing maxillary or mandibular treatment
- Removable and washable cushion cover

### Important

- Can be retrofitted to all units with a comfort or double-jointed headrest
- Double-jointed headrest available in a mechanical, pneumatic or motorised version



## Easy access – thanks to the KaVo armrest retrofit kit for the KaVo ESTETICA E80/E70 Vision and the KaVo ESTETICA E80/E70



With this modern and ergonomically shaped armrest, you will ensure that your patients are in the optimal position and therefore ensure maximum safety during treatment. The right armrest can be folded down completely for easy access. The centred attachment to the patient chair ensures maximum freedom of movement for the practitioner.

### Your benefits:

- More patient and practitioner comfort
- Easier access for elderly and physically disabled patients
- Easy to retrofit

The retrofit kit includes the armrest and the required seat frame and accessories.



## Safety during routine dental practice with the KaVo Original Hoses.

Despite regular care, hoses that are in constant use eventually appear worn out and can make the entire treatment unit look unhygienic. This is why you should ensure that you change them regularly! This way, you guarantee optimal safety for your patients.

Best of all: You receive brand new KaVo hoses at a favourable exchange price!

### Your benefits:

- KaVo original spare parts
- Functionality and quality directly from the manufacturer
- Highest hygiene standards for your dental practice
- Excellent quality at attractive prices
- New look: Increases the attractiveness of your dental practice
- Patient and user safety



## New set of suction hoses for the assistant element


### Your benefits:

- Set of 1x saliva ejector and 1x suction hose (incl. adaptor and two suction cannulas)
- More cost-effective in combination than the individual components
- Fully assembled so that you can replace it yourself
- New material for longer durability and greater functional safety



# Upgrade your treatment unit

Small details can often have a major impact. Increase the treatment comfort for both you and your team and give your patients the good feeling that comes with being in the best hands. Simply upgrade your treatment unit from time to time by upgrading, replacing or repairing parts.

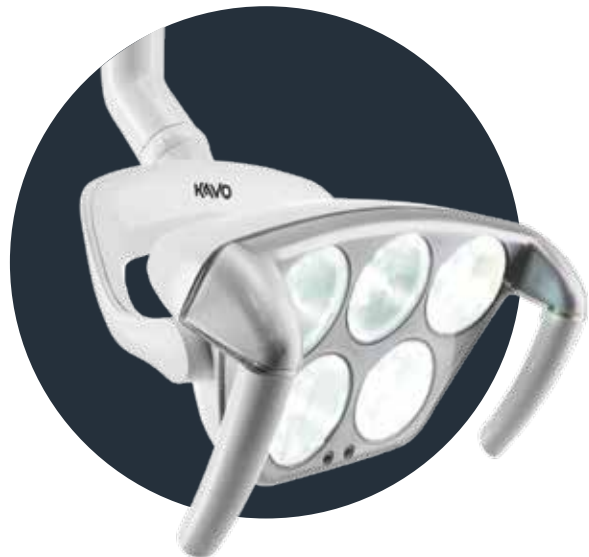


- Light intensity in COMPOsave or dimming mode adjustable from approx. 8,000 to 15,000 lux
- Daylight quality lighting
- Colour temperature adjustable from 4,000 to 6,000 Kelvin
- Silent LED technology
- Removable handles, can be sterilised



Neutral light mode



No-contact light mode change

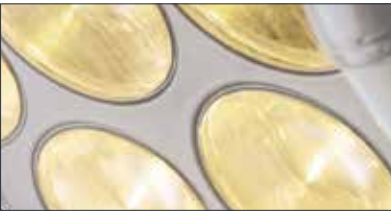

COMPOsave™ mode


## KaVoLUX 540 LED - the perfect light for the highest demands

The KaVoLUX 540 LED can be operated contact-free and therefore hygienically via the sensor or via the intuitive doctor's element. It can also be switched on and off via the automated treatment unit positions. Since it is perfectly suited to your KaVo treatment unit, the KaVoLUX 540 LED offers the perfect hygiene standards, the highest KaVo quality, a long service life and maximum treatment comfort.

A common problem: If the treatment light is not dimmed, the composite often hardens too quickly. This is not the case with the KaVoLUX 540 LED – the innovative COMPOsave™ mode, which is compatible with light-activated restorative materials, filters out all blue light components, slowing down the curing of the filling material. With a brightness of 15,000 lux, you always have the best view for optimal treatment and excellent results.

# All's well that ends well


Three attractive designs can be combined with any upholstery colour:

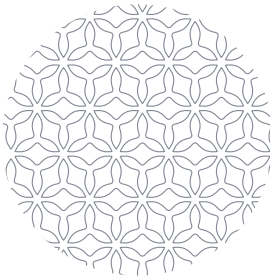

Lisboa Design



Marrakech Design

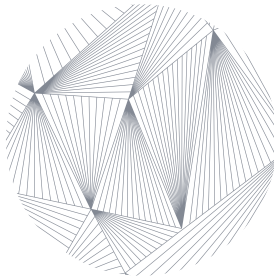

New York Design



## Benefits

- Prolongs your pleasure from your KaVo dental unit
- Increases the attractiveness of your dental practice for patients
- Aesthetic and affordable upholstery protection
- Simple to affix and easy to clean

## Available for the following treatment units:

E80/E80 Vision, E70/E70 Vision, E50/E50 Life, 1058/1058 Life.



## The wireless foot pedal for KaVo ESTETICA E80/E70/E50 makes hygienic treatment a matter of course.

- No unnecessary touching of operating elements, which prevents cross-contamination
- It is possible to operate all instruments safely due to centring – instruments start at a stored speed

- The full speed range is always available
- Centring is not active in level E
- High treatment safety due to indication-related storage of speeds
- All it takes is a light touch on the KaVo foot pedal buttons: Operation of all KaVo ESTETICA treatment unit series instruments, chairs and multimedia functions
- A special control valve for intensity control, e.g. of the turbine compressed air and all other instruments, for smooth and precise instrument control
- Speed control, level switching, centring, spray preselection, blast air, motor direction, chair and multimedia functions can be operated identically
- Best possible adaptation to learned, familiar working methods
- Greater safety due to the absence of a connecting cable between the treatment unit and the foot pedal
- Ergonomic working posture due to free, flexible positioning





### Cannulas for KaVo multifunction and 3-function handpiece

Replacement cannulas for the KaVo multifunction and 3-function handpiece, can be sterilised up to 138 °C (280 °F).

**Your benefits:**

- Always enough replacement cannulas during the running thermal and sterilisation cycles
- Cost-effective replacement of unsightly cannulas
- Easy to remove – simply peel off and replace
- Reusable cannulas are eco-friendly

- For better identification, the cannulas on the multifunction handpiece and on the 3-function handpiece have different labels:
- Golden ring: Cannula for the multifunction handpiece
- Blue ring: Cannula for the 3-function handpiece

Only cannulas affixed with the thermal disinfection and sterilisation symbol can be reprocessed automatically.

## The KaVo cleaning products

### KaVo OXYGENAL™ 6

**Effective water germ reduction for treatment units.**

OXYGENAL 6 consists of hydrogen peroxide, the effects of which are synergistically enhanced by silver ions. In addition to numerous expert reports confirming the effectiveness of OXYGENAL 6, it is particularly well tolerated by materials and does not lead to any wastewater contamination, since the silver is not present at concentrations that would cause toxicity.

At the concentrations used, OXYGENAL 6 does not pose any danger to patients or staff, since the hydrogen peroxide it contains decomposes into water and oxygen without leaving any residues.



### KaVo Patient Service Box More treatment comfort for your patients.



This patient service station can be easily retrofitted to the light support column. It includes a cup holder, a paper dispenser and a shelf. It can be easily disinfected by wiping due to its surface.

**Your benefits:**

- More patient and treatment comfort at your treatment unit
- Fixed storage space for your patients' personal belongings, e.g. glasses, mobile phone, braces box, keys
- The box is modern, hygienic and safe



### KaVo DEKASEPTOL™ Gel

**Disinfectants and cleaning agents for dental suction systems and treatment units.**

DEKASEPTOL Gel completely coats the suction hoses and systems and clings to critical points, instead of just being flushed through. The gel is mixed ready for use and, thanks to the dispenser, is straightforward and easy to apply. This saves you valuable time and also avoids the product coming into contact with skin or work surfaces.

DEKASEPTOL Gel is bactericidal, fungicidal, virucidal, tuberculocidal, effective against HBV and HIV and therefore provides all-round protection for patients and dental practice staff.

# INSTRUMENTS

INDISPENSABLE  
DOWN TO THE  
SMALLEST DETAIL





# KaVo Original Factory Repairs\*


Value preservation, hygiene, sustainability and the Medical Devices Ordinance: There are many aspects that need to be considered when maintaining and repairing medical devices. We support you with KaVo Original Factory Repairs. Ensure optimal treatment quality, legal security, the best possible reliability and increased longevity – all at attractive fixed prices.

**KLAUS ÜBELHÖR**  
*KaVo Technician*

” I know every detail of our products; I have more than 40 years of experience.

Global resources, such as raw materials, water and land, are becoming increasingly scarce. On the other hand, modern consumers are accustomed to simply disposing of damaged equipment instead of having it repaired – the throwaway mentality is now relatively widespread. The earlier we dispose of equipment during its life cycle, the greater the environmental impact and resource consumption in relation to the usefulness of the equipment. The subsequent increased equipment production consumes energy and causes climate-changing greenhouse gas emissions. This is one of the reasons that the careful use of resources is one of the central challenges of our time.

\*This offer is not available in all countries.  
You can find more information at [www.kavobox.com](http://www.kavobox.com).

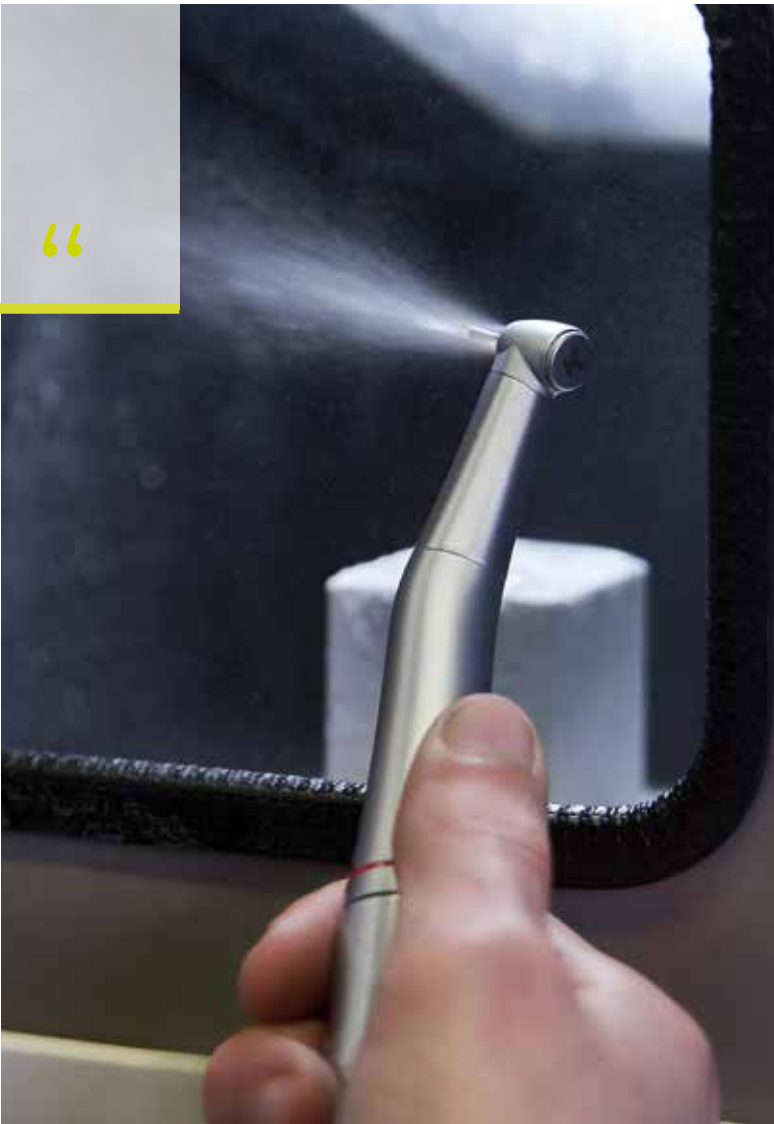

**PETER SCHÖCHLIN**  
*Head of KaVo Original Factory Repairs*

” Your instrument will be almost as good as new again once we have repaired it, and will be guaranteed with a one-year repair warranty. In addition, at least two-thirds of the instrument remains intact during the repair, which conserves resources and is eco-friendly!

### An investment in quality


Anyone who values the highest levels of quality, durability and reliability of instruments, for example, will choose instruments from a premium manufacturer's portfolio, such as KaVo. Anyone who also values maximum value preservation will regularly entrust them to a professional service, such as KaVo Service. This ensures that premium standards are applied even when a repair is due after longer and more intensive use. This repair is either performed directly by designated specialists directly at KaVo or by certified service partners who only use original parts.

### Reliable and hygienic

The dentist Jochen Kania, who has used KaVo products in his practice for many years, agrees. When procuring new instruments, he attaches great importance to durability and hygiene. “The first question I ask myself when making a purchase is whether the instrument or device meets all the hygiene requirements. This is a problem with many instruments and devices from the low-cost sector, as many are not even authorised for use here in Germany. This is why I personally would never buy a no-name product. KaVo offers a quality that cannot be matched by any low-cost products”. He also holds this view in terms of service, particularly

once he has invested money in high-quality products. Anyone who desires quality would be doing themselves a disservice by cutting costs in terms of servicing. He therefore always relies on the KaVo repair service, or appropriately trained service partners. This way, he can be sure that original parts are used whenever anything needs to be replaced. There are also no unpleasant surprises when repairs are carried out at KaVo, because the repairs are conducted at favourable and transparent fixed prices. “My repairs have been carried out



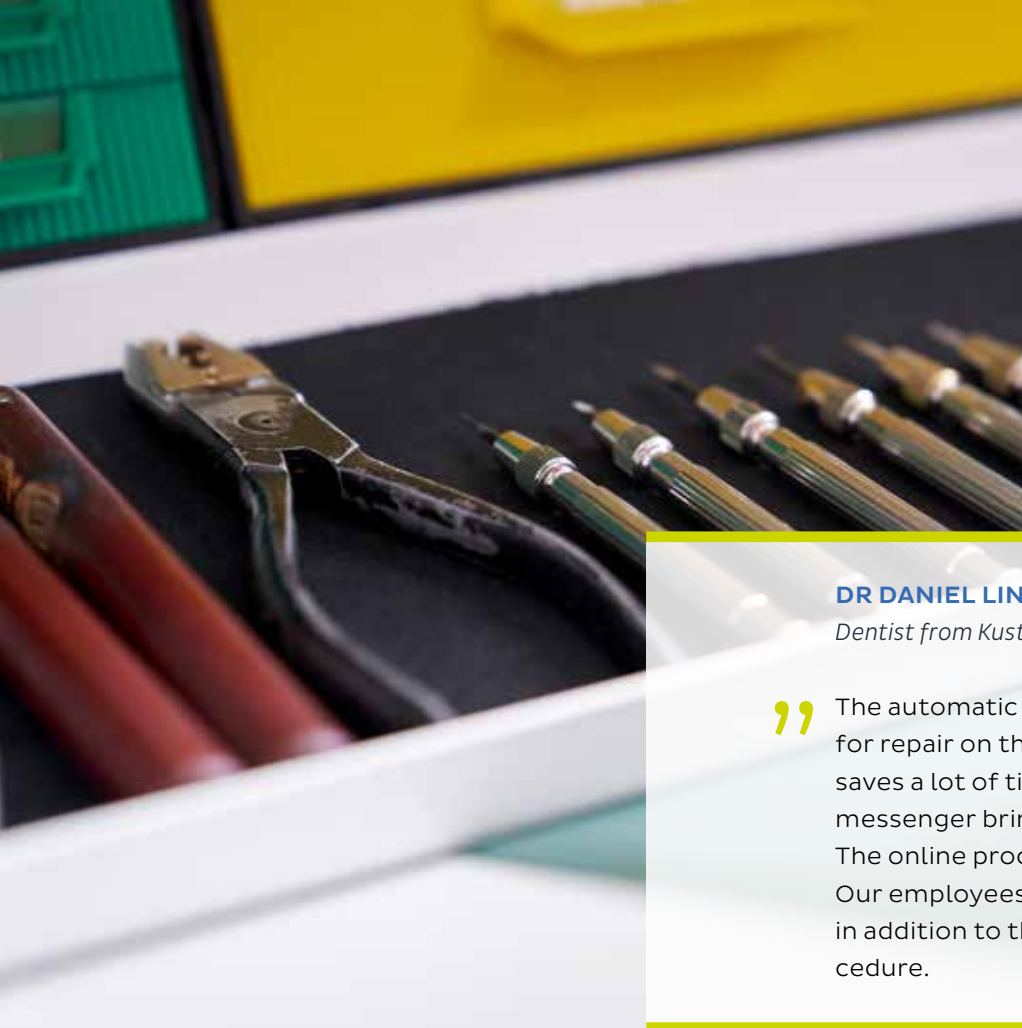

#### DR DANIEL LINDEL

Dentist from Kusterdingen, Germany

” The automatic collection of the parts for repair on the day after the order saves a lot of time, especially as the messenger brings the required sticker. The online procedure is also practical. Our employees can log the repair online in addition to the normal operating procedure. “

exclusively at KaVo for the past one and a half years. A major reason for this – apart from the fixed price guarantee – “is that the process has become much easier and faster with the KaVo BOX, which is why the original factory repair from KaVo is now more interesting for me”, says Kania.

#### Original parts offer safety

In addition to the reliability of the KaVo service, the fact that original parts are exclusively used for repairs plays an important role for the user, as this allows us to guarantee maximum durability, reliable function and, last but not least, absolute legal security. Christian Wiech, who responsible for product management at KaVo, confirms this: “Our service guarantees customers the best service and repair quality for their instruments, since we exclusively use original spare parts and no generic spare parts as a matter of principle”. He also mentions that service becomes a very straightforward and seamless process overall, thanks to the KaVo BOX. “The customer orders our robust KaVo BOX, packs their instruments securely and triggers the repair process. Ordering the KaVo BOX is free of charge, and since it is reusable, we are also contributing to more sustaina-



You can always get it **cheaper**.



You will not get it **better**.



#### DR BERND RIEGER

Dentist from Wiesbaden, Germany

” The upholstered box that I requested free of charge from KaVo was collected on a Monday. I already received the repaired handpiece back on Thursday, and it was clean and tidy. However, my favourite benefit is the fixed price. “

bility through this process. By the way: As long as the product is still under warranty, KaVo will cover the entire cost of the repair”, says Christian Wiech.

#### “Instruments are the interface between the head and the arm”.

Dr Kania cites another reason for his complete reliance on KaVo instruments and KaVo service for regular maintenance: “The direct economic and financial aspect is most important to me. The predefined workflow of the original factory repair also impressed me. After all, my instruments are the interface and the connection between my head and arm. They have to work as one unit. As this is a manual activity, the quality of the instrument is of utmost importance to me. This is why I only use very high-quality products”.

High-quality, and therefore long-lasting, products come at a price – another argument for preserving their value for as long as possible. This is also important to Jochen Kania, who exclusively entrusts his instruments to KaVo specialists or certified partners. “The direct economic and financial consideration are at the forefront for me. With factory repairs, I can hand over the


responsibility because the costs are clear in advance. The process around the KaVo BOX has a pleasant and well-defined workflow and does not require me to perform any approval. The fact that the factory repair process also takes into account sustainability by using limited resources more responsibly is an added bonus”.



You know in advance how much you will save later:  
**Factory repair at attractive and fixed prices\***



[www.kavobox.com](http://www.kavobox.com)

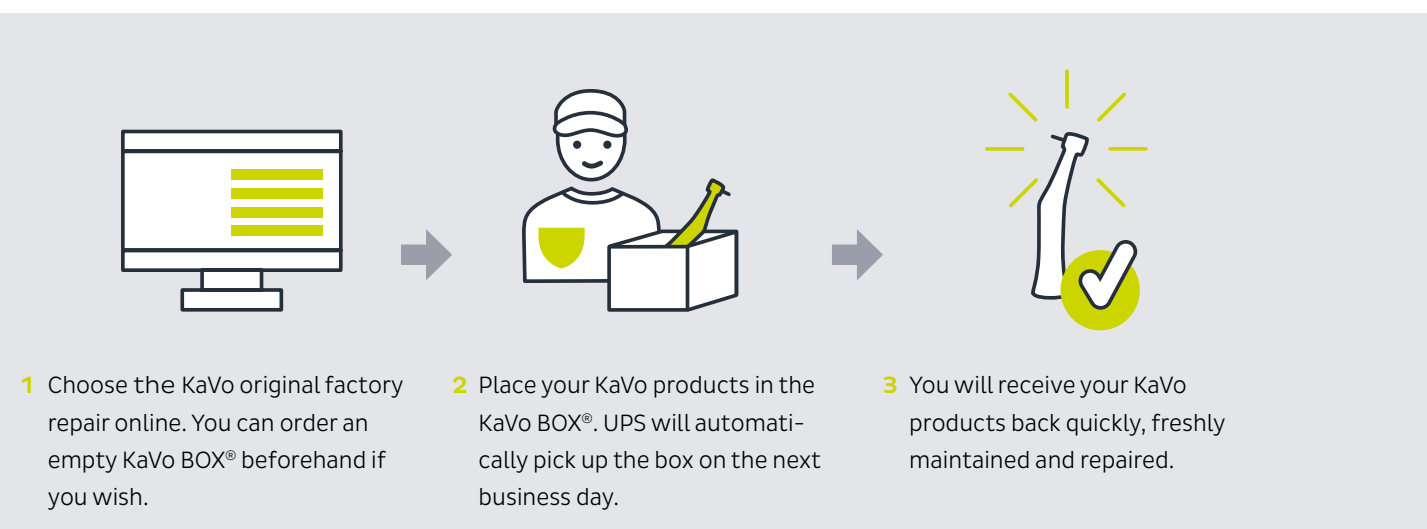

The KaVo BOX® by KaVo offers an attractive and customer-friendly service offer and set standards. KaVo BOX® is the online portal for KaVo original factory repairs.

**Fast, comprehensive and simple repair processing for almost all KaVo instruments and small devices.**

- Maximum cost transparency thanks to attractive fixed prices: You know in advance how much you will save later
- A warm welcome awaits you: Automatic € 30 bonus upon initial registration
- Stay up-to-date: Real-time information on the status of your repair at any time
- The story so far: A history of all repair orders
- At the touch of a button: Free automatic collection of the repair within 1 – 2 working days

**Our best performance for dentists and dental technicians:**

- Safety: Compliance with CE standards
- Value preservation: Technically, the repaired product is as new
- KaVo original factory repair is performed by us directly, without diversions
- Specialised technicians with a high level of experience and competence
- Maximum quality through the use of KaVo original spare parts, tools and testing equipment
- Fast, free collection and return
- Attractive fixed prices for almost all repairs
- Full cost coverage within the warranty period
- 12-month warranty on repairs and replaced parts
- Ordering online saves you time

\*This offer is not available in all countries.  
You can find more information at [www.kavobox.com](http://www.kavobox.com).

**Now available:**  
KaVo BOX® PLUS 2  
Repair cost protection\*



Full repair cost protection for **2 years** after repair.

**Your benefits:**

- 2-year full warranty after completion of repair
- Take a look at the workflow at [www.kavobox.com](http://www.kavobox.com)

**Prerequisites:**

- Your instrument is a maximum of 8 years old
- Only in conjunction with a repair
- Not completely damaged, no maintenance defects, no generic spare parts already installed



Select your product for repair in the repair shop with the KaVo BOX® PLUS symbol.



[www.kavobox.com](http://www.kavobox.com)



\*This offer is not available in all countries.  
You can find more information at [www.kavobox.com](http://www.kavobox.com).



## KaVo replacement motors\*: High power at a low price

As part of our replacement process, you receive KaVo motors that are as good as new directly from us, the manufacturer. Completely overhauled and excellently reconditioned!

### Your benefit:

- Highest quality and long service life
- Cost-effective due to the replacement process
- Quickly available
- Resource-sparing

### Important:

- In addition to thorough cleaning, all wearing parts are replaced with new parts

\*This offer is not available in all countries.  
Contact us at [service.repairs@kavo.com](mailto:service.repairs@kavo.com).



## QUATTROcare PLUS Good care pays off

Your instruments feel like new thanks to care with QUATTROcare PLUS. The validated care process and the ideal amount of oil act like a spa: The interior values are superbly maintained and the instruments are perfectly prepared for routine dental practice.

### KaVo QUATTROcare PLUS 2124 A

Suitable for 3x INTRAMatic short and 1x MULTIflex care coupling, including 1 can of QUATTROcare Plus spray. Care couplings for all common competitor products are available as accessories.

- Just a one-minute care cycle for 4 instruments
- Program-controlled, automatic care cycle
- Optimum oil care thanks to a single-user automatic dosing system
- Easy handling prevents maintenance errors
- Purging of gear and spray channels thanks to a purging function
- No residual oil leakage after care
- Integrated collet chuck care program
- Basic equipment with 4 maintenance couplings (1x MULTIflex/3x INTRAMatic short)
- Universal usability (adaptors available for all common makes)
- Modern and ergonomic design
- 12-month warranty



### KaVo QUATTROcare Plus Spray

Care spray for use in KaVo QUATTROcare PLUS care systems, eco-friendly propellant, each containing 500 ml content



# KaVo spray portfolio for manual processing of your dental instruments



**KaVo CLEANspray/DRYspray**  
Cleaning and drying spray for manual, validated, intensive internal cleaning of instruments, with an expert certificate indicating residual protein removal.



**KaVo Spray**  
Universal care spray for manual care of KaVo instruments, eco-friendly propellant, contains 500 ml.






















**INTRA Instrument Stands**  
The new INTRA instrument stand offers: Improved hygiene thanks to replaceable cell support, better overview of the instruments, practical handling for faster dental practice procedures and improved safety.



KaVo treatment units

Maintenance Kit	
1.010.7974	Maintenance kit, water inlet, compact
1.0108.043	Maintenance kit, water inlet, DVGW
Extended warranties <span>Please take note of the KaVo warranty conditions</span>	
1.010.1621	KaVo PLUS 1 year – extended warranty for E30
1.010.1627	KaVo PLUS 2 years – extended warranty for E30
1.010.1632	KaVo PLUS 4 years – extended warranty for E30
1.010.1624	KaVo PLUS 1 year – extended warranty for 1058/E50/DSE
1.010.1629	KaVo PLUS 2 years – extended warranty for 1058/E50/DSE
1.010.1634	KaVo PLUS 4 years – extended warranty for 1058/E50/DSE
1.010.1626	KaVo PLUS 1 year – extended warranty for E70/E80
1.010.1631	KaVo PLUS 2 years – extended warranty for E70/E80
1.010.1636	KaVo PLUS 4 years – extended warranty for E70/E80
Service Protect+ Packages <span>Please take note of the KaVo warranty conditions</span>	
<b>1.013.5075</b>	<b>KaVo Service Protect+ E70/E80 Vision includes:</b>
1.010.1631	KaVo PLUS 2 years – extended warranty for E70/E80
0.489.3451	Oxygenal 6 (12 x 1L)
1.004.6827	DEKASEPTOL Gel Comfort Set (12 x 1L)
1.005.7679	Footrest protector E70/E80
<b>1.013.5076</b>	<b>KaVo Service Protect+ 1058/E50 Life includes:</b>
1.010.1629	KaVo PLUS 2 years – extended warranty for 1058/E50/DSE
0.489.3451	Oxygenal 6 (12 x 1L)
1.004.6827	DEKASEPTOL gel comfort set (12 x 1L)
1.001.5845	Footrest protector
<b>1.014.6058</b>	<b>KaVo Service Protect+ comp. 1058/E50 Life</b>
1.010.1629	KaVo PLUS 2 years – extended warranty for 1058/E50/DSE
1.010.7974	Maintenance kit, water inlet, compact
1.000.7204	Beginner Set DEKASEPTOL – 2 Pack (3 x 1L)
KaVo Original Suction Hose Set	
<b>0.763.9661</b>	<b>Suction hose set includes:</b>
3.004.7180	Spray mist suction hose, mounted
3.004.7171	Saliva ejector hose, mounted
0.764.5783	Cannula holder, short
0.764.5853	Cannula holder, long
0.764.5863	Cannula holder, small
0.764.5873	Cannula adaptor
0.764.5883	Cannula adaptor
0.236.2129	Suction cannula set

KaVo Original New Upholstery Set			
1.013.8000	All	Progress backrest   Comfort backrest Standard upholstery   RELAXline upholstery   Colours	
Standard Upholstery			
0300801	ESTETICA E70/E80 Vision ESTETICA E70/E80	• Progress backrest • Comfort backrest	
0300803	ESTETICA E50/E30   Primus 1058 Life ESTETICA E50 Life   ESTETICA E50 Primus 1058 Life /Primus 1058 With Standard patient chair	• Double-jointed headrest with rotary control knob • Double-jointed headrest (not available for E30)	
0300807	ESTETICA E50   Primus 1058 Life ESTETICA E50   Primus 1058 With COMPACTchair	• Progress backrest • Comfort backrest • Double-jointed headrest with rotary control knob • Double-jointed headrest	
0300805	ESTETICA Sensus 1066 ESTETICA Comfort 1065	• Somatic headrest (ball pin) • Motor driven or double jointed headrest (socket on cushion) • Somatic headrest (socket on cushion) • Progress backrest • Comfort backrest	
0300805	Status 1080* Globus 1078	• Standard headrest • Somatic headrest (ball pin) • Motor driven or double jointed headrest (socket on cushion) • Somatic headrest (socket on cushion) • Progress backrest • Comfort backrest	
RELAXline Soft Upholstery			
0300802	ESTETICA E70/E80 Vision	• Progress backrest • Comfort backrest	
0300804	ESTETICA E50 Life/ Primus 1058 Life   ESTETICA E30	• Double-jointed headrest with rotary knob • Double-jointed headrest (not available for E30)	
<div><div> Double-jointed headrest with rotary control knob</div><div> Double-jointed headrest</div><div> Standard headrest</div><div> Somatic headrest (ball pin)</div><div> Motor driven or dou- ble-jointed headrests</div><div> Somatic headrest (socket on cushion)</div></div>			
PHYSIO			
0300806	PHYSIO Evo/Evo F/Evo One*		
Standard upholstery colours		RELAXline Soft upholstery colours	
 black no. 33	 pearl grey no. 45	 purple no. 52	 cashmere no. 110
 midnight blue no. 39	 mint no. 38	 berry no. 60	 graphite no. 111
 smoke blue no. 66	 greenery no. 70	 ruby red no. 63	 agave no. 112
 ocean blue no. 64	 emerald green no. 69	 orange no. 59	
 anthracite no. 46	 chocolate brown no. 62		

\* This model may be subject to additional costs due to technical modifications. An additional tubular frame 1.003.9179 and a footrest protector 1.001.5845 may have to be ordered.  
\*\* Also available as RELAXline soft upholstery.



Treatment Unit Hose Overview

ESTETICA E70/E70 Vision | ESTETICA E80/E80 Vision

Intake				
	Doctor's element - S (basic grey)		Doctor's element - T (basic grey)	
	1300 mm	1600 mm	1300 mm	1600 mm
Turbine hose	1.007.8151		1.004.4638	1.004.8367
Motor hose KL701/702/703	1.007.8149		1.004.4635	1.004.8366
Motor hose KL702 ENDO			1.005.1893	
COMFORTbase 404L	1.007.8204		1.004.9811	
Syringe hose MF   3F	1.007.8152 (Doctor)	1.004.8368 (Assistant)	1.004.4654 (Doctor)	1.004.8368 (Assistant)
Piezo-Lux			1.004.4653	
Piezo-LED	1.007.4002		1.007.4002	
Satelec Mini L.E.D.			1.004.8826	1.003.8120 (Assistant)
HF cable			1.004.7137	
ERGOCam4 tubing	1.005.3277	1.005.6302	1.005.3277	1.005.6302
Surgical motor hose SL550			1.001.2651 (2000 mm)	1.004.6825 (3000 mm)
Surgical motor hose SL600			1.009.1700 (2000 mm)	

Suction		
	Length (in mm)	S   T (basic grey)
Spray mist hose, complete	1600	3.004.7180
Spray mist hose, short	1600	1.013.5157
Spray mist hose, roll	25000	1.013.5158
Saliva ejector hose, complete	1600	3.004.7171
Saliva ejector hose, short	1600	1.013.5156
Saliva ejector hose, roll	25000	1.013.5155
Hydrocolloid tube	1500	1.001.2424

Version: KaVo\_HoseOverview\_2021-04-14  
Last updated: 26. October 2020

ESTETICA E50/E50 Life | Primus 1058 Life

Intake				
	Doctor's element - S (basic grey)		Doctor's element - T (basic grey)	
	1300 mm	1600 mm	1300 mm	1600 mm
Turbine hose	1.007.8151		1.004.4638	1.004.8367
Motor hose K200	1.007.8150		1.004.4637	
Motor hose KL702/703	1.007.8149		1.004.4635	1.004.8366
COMFORTbase 404L	1.007.8204		1.004.9811	
Syringe hose MF / 3F (Doctor)	1.007.8152		1.004.4654	
Syringe hose MF / 3F (Assistant)	0.516.8090 (Assistant 1400 mm)	0.516.8190 (Assistant 1600 mm)	0.516.8090 (Assistant 1400 mm)	0.516.8190 (Assistant 1600 mm)
Syringe hose 3F Assistant (r, l)		1.009.0980 (Assistant 2400 mm)		1.009.0980 (Assistant 2400 mm)
Piezo-LED	1.007.4002		1.007.4002	
Satelec Mini L.E.D.		1.003.8120 (Assistant)		1.003.8120 (Assistant)
Satelec Mini L.E.D. Assistant (r, l)		1.009.0979 (Assistant 2400 mm)		1.009.0979 (Assistant 2400 mm)
HF cable			1.004.7137	
ERGOCam4 tubing	1.005.3277	1.005.6302	1.005.3277	1.005.6302

Suction		
	Length (in mm)	S   T (basic grey)
Spray mist hose, complete	1600	3.004.7180
Spray mist hose, complete (r, l)	2400	3.004.9405
Spray mist hose, short	1600	1.013.5157
Spray mist hose, roll	25000	1.013.5158
Saliva ejector hose, complete	1600	3.004.7171
Saliva ejector hose, complete (r, l)	2400	3.004.7136
Saliva ejector hose, short	1600	1.013.5156
Saliva ejector hose, roll	25000	1.013.5155
Hydrocolloid tube	1500	1.001.2424

Version: KaVo\_HoseOverview\_2021-04-14  
Last updated: 26. October 2020

Treatment Unit Hose Overview

1063 | 1065 | 1066 | 1080 | 1058

Intake							
Doctor's element	Abbreviation	1063/65 (basic grey)	1066 (basic grey)		1080 (basic grey)		1058 (basic grey)
		1250 mm	1400 mm	1250 mm	1250 mm	1400 mm	1250 mm
	on tubing			optional		optional	
Turbine hose	LUIN	0.516.5310	0.516.5350	0.516.5310	0.516.5310	0.516.5350	0.516.5310
Turbine hose, flexible		1.002.7075	1.003.2745	1.002.7075	1.002.7075	1.003.2745	1.002.7075
Motor hose 192	KMOL	0.516.5800			0.516.5800		
Motor hose KL 700/701/702/703	KLMO	0.516.5810	0.516.5900	0.516.5810	0.516.5810	0.516.5900	0.516.5810
Motor hose K190/K196/K200	K200	1.000.5652	1.000.5710	1.000.5652	1.000.5652	1.000.5710	1.000.5652
COMFORTbase 405L							1.006.1211
Syringe hose MF	4-F	0.516.8140	0.516.8090	0.516.8140	0.516.8140	0.516.8090	0.516.8140
Syringe hose 3F	4-F	0.516.8140	0.516.8090	0.516.8140	0.516.8140	0.516.8090	0.516.8140
Sonosoft-Lux/ Piezo-Lux	SONO	0.516.1060	0.516.1110	0.516.1060	0.516.1060	0.516.1110	0.516.1060
Piezo-LED							1.007.4002
HF cable		0.226.3832	0.226.3834	0.226.3832			

Suction				
Assistant's element	Abbreviation	Length (in mm)	1063/65   1066   1080 (basic grey)	1058 (basic grey)
	marked			
POLYlux 2	SONO	1550	0.516.1090	
Satelec Mini L.E.D.		1600	1.003.8120	
Hydrocolloid tube		1500	1.001.2424	
Spray mist hose, complete		1600	3.004.7180	
Spray mist hose, complete (r, l)		2400	3.004.9405	
Spray mist hose, short		1600	1.013.5157	
Spray mist hose, roll		25000	1.013.5158	
Saliva ejector hose, complete		1600	3.004.7171	
Saliva ejector hose, complete (r, l)		2400	3.004.7136	
Saliva ejector hose, short		1600	1.013.5156	
Saliva ejector hose, roll		25000	1.013.5155	
Syringe hose 3F	LUIN	1550	0.516.8190	
Syringe hose MF	LUIN	1550	0.516.8190	
Syringe tubing 3F (assistant) (r, l)		2400	1.009.0980	
ASSISTlux hose complete		1600	1.002.1907	
ASSISTlux hose		1600	1.001.4011	
Turbine hose 1058 O	LUIN	1550		0.516.5340
Turbine hose 1058 O (flexible)	LUIN	1550		1.003.2894

Version: KaVo\_HoseOverview\_2021-04-14  
Last updated: 26. October 2020 (Extranet hose numbers overview)



Treatment Unit Hose Overview





ESTETICA E30

Intake				
	Doctor's element - S (basic grey)		Doctor's element - T (basic grey)	
	1300 mm	1600 mm	1300 mm	1600 mm
Turbine hose	1.007.8151		1.007.8151	
Turbine hose, Alu	3.005.5450		3.005.5450	
Motor hose KL701/703	1.007.8149		1.007.8149	
Syringe hose 3F (0.773.0500)	1.007.8152 (Doctor)		1.007.8152 (Doctor)	
Syringe hose 3F (0.773.0500)		0.516.8190 (Assistant   1600 mm)		0.516.8190 (Assistant   1600 mm)
Syringe 3F One incl. Hose	1.009.8741 (Doctor   1300 mm)		1.009.8741 (Doctor   1300 mm)	
Syringe 3F One incl. Hose		1.009.8742 (Assistant   1600 mm)		1.009.8742 (Assistant   1600 mm)
Piezo-LED	1.007.4002		1.007.4002	
Satelec Mini L.E.D. Assistant (r, l)	1.003.8120		1.003.8120	

Suction		
	Length (in mm)	S   T (basic grey)
Wet suction/direct discharge		
Spray mist hose, complete	1600	3.004.7175
Spray mist hose (hose only)	1600	1.013.5157
Spray mist hose, roll	25000	1.013.5158
Saliva ejector hose, complete	1600	3.004.7155
Saliva ejector hose (hose only)	1600	1.013.5156
Saliva ejector hose, roll	25000	1.013.5155
Venturi pump		
Saliva ejector hose, Venturi	1300	1.005.4634
Spray mist tube, Venturi	1300	1.005.4633

Version:  
KaVo\_HoseOverview\_2021-04-14

Last updated: 26. October 2020

Footrest protector		
Footrest protector neutral		
1.001.5845	Footrest protector 1058/E50	
1.005.7679	Footrest protector E70/E80	
1.013.6513	Footrest protector RELAXline 1058/E50	
1.013.6514	Footrest protector RELAXline E70/E80	
Footrest protector Design		
1.013.0161	Footrest protector E70/E80 New York	
1.013.0164	Footrest protector E30/1058/E50 New York	
1.013.0162	Footrest protector E70/E80 Marrakech	
1.013.0165	Footrest protector E30/1058/E50 Marrakech	
1.013.0163	Footrest protector E70/E80 Lisboa	
1.013.0166	Footrest protector E30/1058/E50 Lisboa	
Retrofit kit, Vision armrest for E70 and E80 Classic		
1.011.5785	Retrofit kit, armrest, left	
1.011.5785	Retrofit kit, armrest, left and right	
Retrofit kit, comfort headrest		
1.000.6654	Comfort head cushion	
Suitable for:	The double-jointed headrest or motorised headrest. (rotating cushion socket on upholstery 1.004.9162)	
	Somatic rotating cushion (sockets on upholstery 0.790.4841)	
	Somatic rotating cushion (ball pin on upholstery 0.790.3651)	
	<b>Attention!</b> Only possible in combination with plate 1.001.4343!	
KaVo foot pedal/wireless foot pedal		
1.006.1271	Wireless foot pedal for E80/E70 Classic	
1.011.0414	Wireless foot pedal for E80/E70 Vision	
1.008.0370	Wireless foot pedal for E50 Life/Classic	

Patient Service Box	
1.002.1067	Patient Service Box
Hygiene products	
0.224.6876	Outlet tap
0.489.3451	Oxygenal 6 1000 ml
1.000.7204	Beginner Set Dekaseptol 3 x 1 L
1.000.7205	Refill canister Dekaseptol 6 L
1.000.7206	Dosing Dispenser Dekaseptol
1.004.6827	DEKASEPTOL gel comfort set (12 x 1L)
1.007.2427	Hygiene package
KaVo 540 LED	
1.009.4271	KaVoLUX 540 LED T R (for right-hand unit)
1.009.4272	KaVoLUX 540 LED T L (for left-hand unit)
1.009.4273	KaVoLUX 540 LED T R/L (for right/left swivel units)
1.009.1430	KaVoLUX 540 LED T R Centro (for mounting on Centro)
1.008.7973	KaVoLUX 540 LED T R ceiling mounting
1.008.7974	KaVoLUX 540 LED T R/L ceiling mounting (for right/left swivel units)
1.008.8323	KaVoLUX 540 LED U R Centro (mounting on Centro)
1.009.4266	KaVoLUX 540 LED U R (for right-hand unit)



# KaVo Instruments

Replacement equipment*	
0.999.0178	A-INTRA K-LUX motor 196 0676.0800
0.999.0179	A-Intra-K motor 197 0676.0810
0.999.0182	A-valve body, mounted 0773.0021
0.999.0413	A-K-LUX motor 192 0676.0650
0.999.0457	A-SONOsoft LUX 0571.5421
0.999.0550	A-INTRAmatic LUX motor 700 KL 0676.1000
1.000.4137	A-PIEZOlux handpiece 0571.5211
1.005.4943	A-INTRA K-LUX motor K200 1000.6717
1.005.4945	A-INTRA K-LUX motor K201 1001.8156
1.005.4947	A-INTRA LUX MOTOR KL 701 1001.8776
1.005.4948	A-INTRAmatic LUX 702KL 1003.5622
1.006.6650	A-MFH valve body, complete 1002.1650
1.008.7015	A-surgical motor INTRA SL 550 1001.3421
1.008.7017	A-surgical motor INTRA S550 1000.8072
1.008.9078	A-INTRA LUX KL 703 LED 1007.0150
1.009.2580	A-valve body 3F-syringe 0773.0261
1.014.8141	A-INTRA LUX S600 LED 1008.8000
KaVo BOX® – Extended warranty**	
Please take note of the KaVo warranty conditions	
KaVo PLUS 1 year**	
1.014.3872	KaVo PLUS 1 year INTRA head
1.014.3873	KaVo PLUS 1 year Coupling
1.014.3874	KaVo PLUS 1 year Contra-angle handpiece gr/bl w/o head
1.014.3875	KaVo PLUS 1 year Handpiece
1.014.3876	KaVo PLUS 1 year Scaler/ZEG
1.014.3877	KaVo PLUS 1 year PROPHYflex/RONDOflex
1.014.3878	KaVo PLUS 1 year Contra-angle handpiece gr/bl w. head
1.014.3879	KaVo PLUS 1 year Turbine
1.014.3880	KaVo PLUS 1 year Contra-angle handpiece red
1.014.3881	KaVo PLUS 1 year Motor
1.014.3882	KaVo PLUS 1 year MF syringe
1.014.3883	KaVo PLUS 1 year 3F syringe
1.014.3884	KaVo PLUS 1 year Laboratory control unit
1.014.3885	KaVo PLUS 1 year Laboratory turbine/handpiece

\* This offer is not available in all countries. Please contact us at [service.repairs@kavo.com](mailto:service.repairs@kavo.com).  
\*\* This offer is not available in all countries. You can find information at [www.kavobox.com](http://www.kavobox.com).



KaVo PLUS 2 years**	
1.014.3886	KaVo PLUS 2 years INTRA head
1.014.3887	KaVo PLUS 2 years Coupling
1.014.3888	KaVo PLUS 2 years Contra-angle handpiece gr/bl w/o head
1.014.3889	KaVo PLUS 2 years Handpiece
1.014.3890	KaVo PLUS 2 years Scaler/plaque and calculus remover
1.014.3891	KaVo PLUS 2 years PROPHYflex/RONDOflex
1.014.3892	KaVo PLUS 2 years Contra-angle handpiece gr/bl w. head
1.014.3893	KaVo PLUS 2 years Turbine
1.014.3894	KaVo PLUS 2 years Contra-angle handpiece red
1.014.3895	KaVo PLUS 2 years Motor
1.014.3896	KaVo PLUS 2 years MF syringe
1.014.3897	KaVo PLUS 2 years 3F syringe
1.014.3898	KaVo PLUS 2 years Laboratory control unit
1.014.3899	KaVo PLUS 2 years Laboratory turbine/handpiece
KaVo PLUS 3 years**	
1.014.3900	KaVo PLUS 3 years INTRA head
1.014.3901	KaVo PLUS 3 years Coupling
1.014.3902	KaVo PLUS 3 years Contra-angle handpiece gr/bl w/o head
1.014.3903	KaVo PLUS 3 years Handpiece
1.014.3904	KaVo PLUS 3 years Scaler/plaque and calculus remover
1.014.3905	KaVo PLUS 3 years PROPHYflex/RONDOflex
1.014.3906	KaVo PLUS 3 years Contra-angle handpiece gr/bl w. head
1.014.3907	KaVo PLUS 3 years Turbine
1.014.3908	KaVo PLUS 3 years Contra-angle handpiece red
1.014.3909	KaVo PLUS 3 years Motor
1.014.3910	KaVo PLUS 3 years MF syringe
1.014.3911	KaVo PLUS 3 years 3F syringe
1.014.3912	KaVo PLUS 3 years Laboratory control unit
1.014.3913	KaVo PLUS 3 years Laboratory turbine/handpiece
KaVo BOX® – Repair cost protection**	
Please take note the KaVo conditions for repair cost protection	
1.014.3842	KaVo BOX PLUS 2 years Head
1.014.3843	KaVo BOX PLUS 2 years Coupling
1.014.3844	KaVo BOX PLUS 2 years Contra-angle handpiece gr/bl w/o head
1.014.3845	KaVo BOX PLUS 2 years Handpiece
1.014.3846	KaVo BOX PLUS 2 years Scaler/plaque and calculus remover
1.014.3847	KaVo BOX PLUS 2 years PROPHY-/RONDOflex

\*\* This offer is not available in all countries. You can find information at [www.kavobox.com](http://www.kavobox.com).

1.014.3848	KaVo BOX PLUS 2 years Contra-angle handpiece gr/bl w/o head
1.014.3849	KaVo BOX PLUS 2 years Turbine
1.014.3850	KaVo BOX PLUS 2 years Contra-angle handpiece red
1.014.3851	KaVo BOX PLUS 2 years Motor
1.014.3852	KaVo BOX PLUS 2 years MF syringe with cannula
1.014.3853	KaVo BOX PLUS 2 years 3F syringe with cannula
1.014.3854	KaVo BOX PLUS 2 years Laboratory control unit
1.014.3855	KaVo BOX PLUS 2 years Laboratory turbine/handpiece
KaVo BOX® – Fixed price**	
1.013.1019	Small Standard Service Package
1.013.5497	Standard Repair Package HEAD
1.013.5554	Standard Repair Package COUPLING
1.013.5555	Standard Repair Package CONTRA-ANGLE HANDPIECE
1.013.5556	Standard Repair Package HANDPIECE
1.013.5557	Standard Repair package SCALER
1.013.5558	Standard Repair package POWDER JET
1.013.5559	Standard Repair package TURBINE
1.013.5560	Standard Repair package CONTRA-ANGLE HANDPIECE RED
1.013.5561	Standard Repair package MOTOR
1.013.5562	Standard Repair package LABORATORY
1.013.5563	Standard Repair package CONTRA-ANGLE HANDPIECE 2
1.013.5564	Standard Repair package MF SYRINGE
1.013.8504	Standard Repair package 3F SYRINGE
1.013.8120	Standard Repair package LAB CONTROL UNIT
KaVo BOX® – Service Check**	
1.001.8903	Service Check Surgical Unit
1.011.4824	Service Check Diagnosis Handpiece
1.011.5108	Service Check Quattro Care
Cannulas	
0.773.0011	Cannula MF handpiece
0.773.0231	Cannula 3-F-Handpiece
Instrument stands	
3.005.5204	INTRA instrument stands

\*\* This offer is not available in all countries. You can find information at [www.kavobox.com](http://www.kavobox.com).

KaVo CLEANspray   DRYspray	
1.007.0579	KaVo CLEANspray 2110 P package
1.007.0580	KaVo DRYspray 2117 P package
1.007.0573	KaVo CLEAN-/DRYspray starter set 2116 P
1.007.1776	INTRAmatic adaptor
1.007.1775	MULTIflex adaptor
1.007.3350	CLEAN/DRYspray   COMFORTdrive adaptor
0.411.9691	Cleanpac package
KaVo Spray	
0.411.9640	KaVo SPRAY 2112A
1.011.5721	
0.411.9911	INTRA spray head
0.411.9921	MULTIflex spray head
1.005.3154	COMFORTdrive spray head
1.005.8436	NSK spray head, mounted
1.005.8365	SIRONA spray head, mounted

## KaVo e-Service\*

1.013.3410	Installation of CONEXIO software (single workstation)
1.013.3412	Installation of CONEXIO software, 1x Server, 1x Client
1.013.3413	Installation per additional CONEXIO client
1.013.3414	Installation of DIAGNOcam software (single user)
1.013.3415	Installation per additional DIAGNOcam client
1.013.3416	Database relocation to new server
1.013.3408	Connection of 2D software to PMS
1.013.3418	Update within the main version (server only)
1.013.3404	General remote maintenance (per hour)

\*This offer is not available in all countries. Please contact us at [service.multimedia@kavo.com](mailto:service.multimedia@kavo.com).





# Overview of contact partners

## Technical Service

### Instruments

- Handpieces and contra-angle handpieces, turbines
- QUATTROcare
- DIAGNOdent
- MASTERSurg, EXPERTSurg, INTRASurg
- STERlclave
- SPRAYrotor
- HealOzone
- Surgical motor/air motor

Service.Instruments@kavo.com

### Dental technology

- Laboratory drives
- Laboratory devices, laboratory equipment
- Suction systems
- Articulators

Service.Zahntechnik@kavo.com

### CAD/CAM

- ARCTICA
- Everest

Service.cadcam@kavo.com

### CAD/CAM

#### Application Technology

- ARCTICA
- Everest

AWT@kavo.com

### Multimedia

- CONEXIO
- DIAGNOcam
- ERGOcom
- ERGOcam

Service.Multimedia@kavo.com

### Facilities

- Treatment units
- Treatment lights
- Centro
- DOMINO
- PHYSIO
- K/KL Motor

Service.Einrichtungen@kavo.com

### Laser/dental microscopes

- KEY Laser/GENTLEray
- Master lase/Expert lase
- Leica M320

Service.Laser@kavo.com

## Customer Support Center (CSC)

### Repairs

#### Receipt/dispatch

- Orders
- Delivery notes
- Loaned goods

Tel: +49 7351 56-1900  
service.repairs@kavo.com

#### Repairs Processing

- Cost estimates

Tel: +49 7351 56-4900  
service.repairs@kavo.com

**Please always provide the following data when making enquiries:**

- Name of the specialised dealer/technician
- To which type of unit does the enquiry refer (e.g.: ESTETICA E80 Vision)
- What is the serial number of the unit?
- Who is the operator of the unit
  - Practice name
  - Street/house number
  - Postcode/city
- Detailed description of the request

*If the line is busy, please leave a message on our answering machine.*

**We will be happy to call you back.**

## Customer service (order processing)

**Please send your orders by email whenever possible.**

### Contact

#### European Countries

Tel: +49735156-1200  
customerservice-eu@kavo.com

#### Other Countries

Tel: +49735156-1400  
customerservice-global@kavo.com

#### Germany, Switzerland, Austria

Tel: +49735156-1300  
kundenservice-dach@kavo.com



**You are also welcome to use our contact form for enquiries. You can find this and the telephone numbers in the respective countries on our homepage:**

**www.kavo.com**

## Notes

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.



# Dental excellence in every area.



## Practice equipment

KaVo treatment units, treatment lights, dental chairs, patient communication systems and other dental practice accessories.



## Instruments

Dental handpieces and contra-angle handpieces, turbines, powder jet handpieces and small devices for all applications ranging from diagnostics to prophylaxis, and including restoration, surgery, endodontics and instrument care and maintenance.

The products, equipment and services shown and described in this catalogue are not available in all countries. All information corresponds to the knowledge available at the time of printing. KaVo Dental GmbH assumes no responsibility whatsoever for discrepancies in the colour and shape of the images, or for errors or misprints, and reserves the right to make changes to the brochures at any time. Reprinting, even in part, is only authorised with the approval of KaVo Dental GmbH.